

• **Q2 ION** Ultra-Compact Spark-OES Metals Analyzer


Bruker Elemental's all new spark spectrometer Q2 ION elevates metal analysis into new levels of simplicity and ease-of-use. Today Q2 ION is the smallest and lightest ultra-compact spark emission spectrometer for metals analysis available. It is a versatile multi-matrix system for comprehensive incoming material inspection and quality assurance of metal alloys. Its affordable price and low operational costs make it the ideal tool for smaller and medium-size businesses.

Q2 ION covers all major alloying elements in many applications such as Ferrous alloys, Aluminium, Coppers, and many more. It ideally addresses the requirements of small and mid-size foundries, metal processing industries, fabricators, quality control departments, warehouses, metal recyclers, and even inspection companies.



Q2 ION - Metals Analysis made easy

Its design makes Q2 ION ultra light (less than 44lbs/20 kg) so it can easily be hand-carried even to a nearby site for analysis. An optional case is also available. Despite its low weight, it is suitable for applications in which ruggedness is required. Q2 ION also defines new standards in ease-of-use. Place your sample onto the spark stand and press the start button. In less than thirty seconds you get the complete elemental composition of your metal.

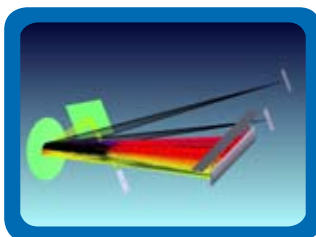


Q2 ION - Patented Optical System

The new patented Flat Field CCD optics is a masterpiece of optics design and mechanical engineering. Active Ambient Compensation (AAC) provides maximum stability in a temperature range between 10 and 45° C (50 and 113° F). The high-definition CCD detector together with well-proven ClearSpectrum® technology provide best-in-class analytical performance.

Detector with ClearSpectrumTechnology



Typical Analysis Screen

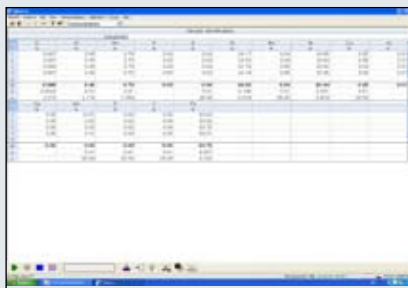

Q2 ION



Technical Data Q2 ION

Patented Optical System	<ul style="list-style-type: none"> Un-coated CCD detector with lowest dark current Flat field grating Full spectrum coverage: 170- 685 nm Resolution: 30 pm Argon purged for best transparency ClearSpectrum® technology for advanced spectra deconvolution Active Ambient Compensation (AAC) for operation between 10 and 45°C (50 and 113° F)
Analytical Solution Packages (ASP)	<ul style="list-style-type: none"> Different matrix calibration packages available ASPs cover all major elements & alloy groups Upgradable for future expansion
Source Generator	<ul style="list-style-type: none"> Maintenance-free, two phase PWM Generator Frequency 50 to 1000 Hz Spark & arc-like discharges from 10 µs to 2 ms
Sparkstand	<ul style="list-style-type: none"> Maintenance-free Argon consumption 2.5 l/min. during measurement Argon quality 4.8 (or better)
Software	<ul style="list-style-type: none"> Intuitive Windows® based software for simple routine operation Various user levels for secure and task-specific operations Functions for qualitative and quantitative analysis Complete QMatrix Software Suite including analysis database & interfaces to Office software Grade Library functions
Electrical Data	<ul style="list-style-type: none"> 100 to 240 V (50/60 Hz) 200 W during measurement, 50 W standby 16 A (240 V) slow blow fuse or 25 A (100 V) slow blow fuse
Dimensions & Weight	<ul style="list-style-type: none"> Width 440 mm / 17 inches Height 220 mm / 9 inches Depth 390 mm / 15 inches Weight ~ 19 kg / ~ 42 lbs

● Bruker Elemental GmbH

Kalkar · Germany
 Phone +49 (28 24) 9 76 50-0
 Fax +49 (28 24) 9 76 50-10
 info@bruker-elemental.com

www.bruker.com/oes

Bruker AXS Inc.

Billerica, MA · USA
 Phone +1 (978) 663-3660
 Fax +1 (978) 667-5993
 elemental@bruker-axs.com